

# Silicone Silo Bag – Clear Cast




## *For the staged reduction of gastroschisis and omphalocele*

Silo bags from SSP provide a secure, closed environment for exposed viscera during the staged closure of congenital ventral wall defects. They help prevent the loss of peritoneal fluid and minimize the risk of bacterial contamination.

Made of clear implantable-grade silicone, our seamless silo bags enable excellent visualization of their contents. Options include wire-spring or wire-free bag openings and a wide-body style (see back), and all bags feature a reinforced silicone loop for continuous traction.

SSP silo bags are indicated for staged bedside reduction or temporary protection prior to O.R. closure.

- ▼ **Seamless, clear cast silicone** elastomer shell allows for better content visualization
- ▼ **WIRE-SPRING or WIRE-FREE** bag opening to be placed in the defect for a closed environment
- ▼ **Wide range of sizes** (see back)
- ▼ **Continuous postoperative traction** with durable silicone loop on distal end
- ▼ **Sizer ring** included to verify desired bag diameter prior to use
- ▼ **Size imprinted** on each bag
- ▼ **Not made** with natural rubber latex
- ▼ **Supplied sterile**

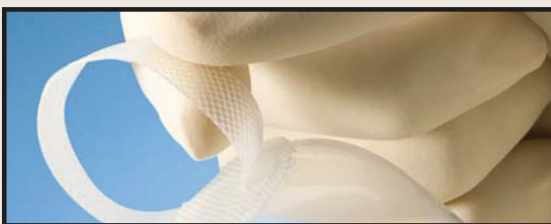

Choose from bag openings with a **wire spring** encapsulated in silicone or a **wire-free** ring. The opening is placed over the organs, gently compressed to permit insertion into the defect and released. It may then be sutured into place. All bags feature a reinforced distal loop for continuous overhead traction.



SPECIALTY SURGICAL PRODUCTS, INC.

# Silicone Silo Bag – Clear Cast

Catalog No.		Description	Diameter	Length
Wire Spring	Wire-Free			
SSB-WR03	N/A	Silicone Silo Bag	3.0cm	8.0 in.
SSB-WR04	N/A	Silicone Silo Bag	4.0cm	8.0 in.
N/A	SSB40	Silicone Silo Bag	4.0cm	8.0 in.
N/A	SSB40W	Silicone Silo Bag, Wide Body	Ring: 4.0cm Silo: 5.5cm	8.0 in.
SSB-WR05	N/A	Silicone Silo Bag	5.0cm	8.0 in.
N/A	SSB55	Silicone Silo Bag	5.5cm	8.0 in.
SSB-WR06	N/A	Silicone Silo Bag	6.0cm	8.0 in.
N/A	SSB65	Silicone Silo Bag	6.5cm	8.0 in.
SSB-WR07	N/A	Silicone Silo Bag	7.0cm	8.0 in.
N/A	SSB75	Silicone Silo Bag	7.5cm	8.0 in.
SSB-WR08	N/A	Silicone Silo Bag	8.0cm	8.0 in.
SSB-WR09	N/A	Silicone Silo Bag	9.0cm	8.0 in.
N/A	SSB95	Silicone Silo Bag	9.5cm	8.0 in.
SSB-WR10	N/A	Silicone Silo Bag	10.0cm	8.0 in.



SSP also offers a wide-body silo bag with a 5.5cm diameter that tapers to a 4.0cm wire-free ring. This provides more room for the contents without enlarging the defect.

**CONTACT CUSTOMER SERVICE FOR CURRENT PRICING**



SPECIALTY SURGICAL PRODUCTS, INC.

1123 NORTH US HIGHWAY 93 + VICTOR, MT 59875 + 888.878.0811 + 406.961.0102 + 406.961.0103 FAX  
WWW.SSP-INC.COM + CUSTOMERSERVICE@SSP-INC.COM

Order for Doctor: \_\_\_\_\_

Patient: \_\_\_\_\_

Surgery Date: \_\_\_\_\_

Catalog No: \_\_\_\_\_

Quantity: \_\_\_\_\_

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

**▼ LIMITED WARRANTY, LIMITATION OF LIABILITY AND DISCLAIMER OF OTHER WARRANTIES**

SSP warrants that reasonable care was used in the manufacture and production of this product. Because SSP has no control over the conditions of use, patient selection, surgical procedure, post-surgical stresses, or handling of the device after it leaves our possession, SSP does not warrant either a good effect or against an ill effect following its use. SSP shall not be responsible for any incidental or consequential loss, damage or expenses directly or indirectly arising from the use of this product. This warranty is in lieu of and excludes all other warranties not expressly set forth herein, whether express or implied by operation of law or otherwise, including but not limited to any implied warranties of merchantability or fitness for use.